English



Specialized Hearing Solutions Model SBJ525SS

Sonic Bomb Alarm Clock & Bed Shaker



OPERATION INSTRUCTIONS

IMPORTANT - Please read these instructions carefully before use and retain for future reference. See operating instructions on back.

- Use only the supplied adaptor to connect the unit,


Warning: To prevent fire or shock hazard do not expose this appliance to rain or moisture. Before connecting to mains make sure that the requirements specified on the apparatus corresponds to the available power supply i.e. AC120V 60Hz.

- -- The ventilation should not be impeded by covering the ventilation openings with items such as, newspaper, table cloth, curtains, etc.
- -- No naked flame sources, such as lighted candles should be placed on the apparatus.
- -- It shall not be exposed to dripping or splashing and that no objects filled with liquids, such as vases, shall be placed on the apparatus.
- -- Do not place the product in closed bookcases or racks without proper ventilation.
- -- Where the mains plug is used as the disconnect device, the disconnect device shall remain readily operable.

GENERAL DATA

Power Consumption: Adaptor AC6V~ 600mA) UL listed transformer (Underwriters Laboratories) for USA

INSTRUCTIONS LED ALARM CLOCK WITH BED SHAKER MODEL SBJ525SS

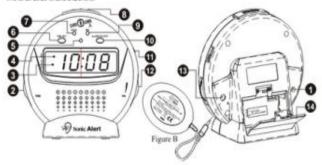

Features show on diagram above:

6. Time Set 1. AC adpator input 11. Function Switch 2. Tone 7. Hour 12. Volume Knob 8. Snooze 3. Alarm indicator 13. Vibrator out Jack

4. PM indicator 9. Minute 10. Alarm Set / Off

14. Battery Compartment door

Setting the correct time -

5. Flashing LED

1) Press and hold the Time Set Button (6) at the same time press the Hour Button (7) to advance the clock display rapidly until the hour is correct.

2) Press and hold the Time Set Button (6) at the same time press the Minute Button (9) to advance the clock display rapidly until the minute is correct. A lighted dot (4) in the upper left of the clock display indicated the time is PM. For AM time indication, the dot does not light up.

Setting the alarm time -

1) Press and hold the Alarm Set Button (10) at the same time press the Hour Button (7) to advance the alarm display rapidly until the hour is correct.

2) Press and hold the Alarm Set Button (10) at the same time press the Minute Button (9) to advance the alarm display rapidly until the minute is correct.

3) Make sure the alarm time is correctly set for AM or PM, the lighted dot (4) indicated PM. Please be sure to release the HOUR and MIN set buttons before releasing the ALARM button.

Waking to Alarm -

- 1) Set Function Switch(11) to choose "VIB", "BUZZ" or "VIB / BUZZ"
- 2) Set desired volume using the Volume knob (12) for buzzsound.
- 3) Set alarm wake time per instructions for setting alarm time above.

Battery Backup -

To use the battery backup, you can install a 9-volt alkaline battery in the back of the clock located under the battery cover (14). Battery not included. Note the LED display will not show during power failure but the clock will still be running and keeping proper time until the power is restored.

NOTE: For the best results it is recommended that the battery be changed every 6-8 months. If the clock is not used for a while, it is recommended to remove the battery to aviod leakage.

Shutting Off Alarm -

Alarm can be shut off temporarily for the day or shut off permantently.

- 1) When the alarm goes off, simply press the AL off button (10) the alarm will shut itself off for the day and come back on the next day at the same time if the alarm setting has not changed.
- 2) To shut off the alarm permanently, slide function switch (11) to off positon. Alarm will not come on until Function Switch (11) is set to "VIB" or "BUZZ" or "VIB/BUZZ" position.

Snooze Operation -

SNOOZE button (8)-press after alarm has sounded. The alarm will be turned off for 9 minutes for extra sleep. The alarm will come back again after 9 minutes. You may repeat many time up to 119 minutes .

Alarm -

Set the alarm time you wish and slide the Function Switch (11) switch to the desired position. Set to VIB for bed shaker, and BUZZ for sound, or VIB/ BUZZ for both.

Adjust the tone by rotating the Tone wheel (2) until you have reach the desired level.

Volume -

To adjust the volume rotate the Volume wheel (12) until you have reach the desired level.

Bed Shaker / Vibrator - Figure B

- 1) Attach the bed shaker jack to the back of the clock at VIBRATOR input (13).
- 2) To turn vibrator slide function switch (11) to VIB or VIB / BUZZ position.
- 3) Place vibrator under mattress or pillow.

AC 6V Power Input -

Attach the AC power supply input to the back of clock labeled AC6V-IN (1). Important: do not confuse AC input with vibrator input



The Lightning flash and arrow head within the triangle is a warning sign alerting you of dangerous voltage inside the product.

CAUTION:

TO REDUCE THE RISK OF ELECTRICAL SHOCK, DO NOT REMOVE THE COVER (OR BACK) FROM THE PRODUCT. NO USER SERVICEABLE PARTS INSIDE. REFER SERVICING TO YOUR SONIC ALERT DEALER



The exclamation mark within the triangle is a warning sign alerting you of important instructions accompanying the product.

WARNING: TO PREVENT FIRE OR SHOCK HAZARD, DO NOT EXPOSE THIS APPLIANCE TO RAIN OR MOISTURE.

- Attention should be drawn to the environmentak aspects of battery disposal.
- this indicating separate colletion for electrical ans electroinc equipment.

SONIC ALERT'S 1 YEAR LIMITED WARRANTY

The Sonic Boom Alarm Clock Model SBJ525ss is warranted against manufacturing defects in materials and workmanship for (1) year from the date of purchase. Within this period Sonic Alert will repair or replace at our option the SBJ525ss without charge for parts and labor. Simply send the SBJ525ss (postpaid) and a copy of your sales slip as proof of purchase to:

Sonic Alert Inc., 1050 E. Maple Road, Troy MI 48083, USA.

Sonic Alert

Made in China Sonic Alert 1050 East Maple Rd. Trov. MI 48083 www.SonicAlert.com 1-888-864-2446 V/TTY

