

PERSONAL AMPLIFIER SPECIFICATION DATA


1-on-1 Conversations • TV Listening • Listening in the Car • Small Groups and much more


Pocketalker® 2.0 Personal Amplifier




The Pocketalker® 2.0 amplifies sound clearly and easily for better understanding.

Personal Amplifier Model PKT E1:

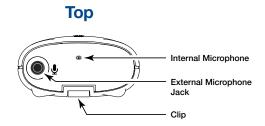

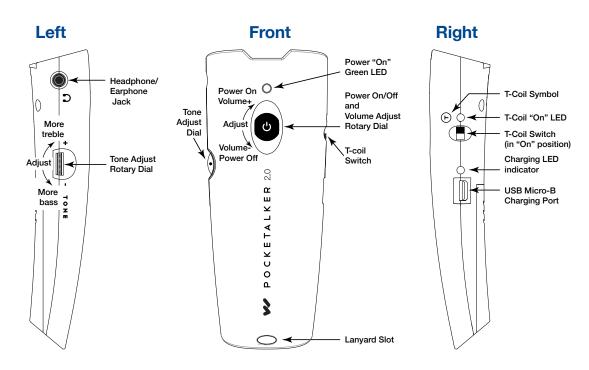

Back

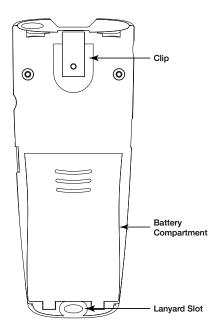

Pocketalker® 2.0 Personal Amplifier

Personal Amplifier Model PKT E1:

Dimensions:	4.9" L x 2.0" W x 1.0" H. (124mm x 51mm x 25mm) (Not including clip or microphone assembly)			
Weight:	•			
Color/ material:	White ABS/Polycarbonate molded plastic case.			
Battery Type:	Two 1.5V (AAA size), Alkaline or NiMH rechargeable*			
	(Williams Sound BAT 010-2 Alkaline or BAT 022-2 NiMH recommended)			
Battery Life:	Up to 105 hours (typical usage) with the BAT 010-2. Low Battery Indicator: Power "On" Green light on front panel flashes.			
**NiMH battery Charging:	blinks when charging.			
Microphone:	Internal, omnidirectional MEMS microphone (micro-electro-mechanical system).			
	Optional External , omnidirectional, plug-in electret microphone with foam windscreen and 3.5mm mono plug (MIC 014-R).			
T-Coil:	Internal with on/off switch and LED indicator.			
Tone Control:	External rotary control knob allows user to adjust high/low frequency to optimize hearing.			
***Balance Control:	Adjusts Left/Right Volume Balance			
Power/ Volume Control:				
Microphone Jack:	· -			
	3.5 mm, stereo mini plug. Tip=Left/Mono, Ring=Right, Sleeve=Gnd.			
Optional Accessories:	Clip: CLP 014; Allows the user to clip MIC 014 to clothing such as a lapel or shirt pocket. Lanyard: RCS 003; Allows user to wear unit around the neck, or as a wrist strap. 3.5mm-3.5mm Extension Cord: WCA 007 WC; extends an external microphone to 12 ft. Stereo Earphones / Headphones / Neckloop: EAR 042 / HED 024, HED 036 / NKL 001-S (stereo) Neckloop Mono Earphones / Headphones / Neckloop***: EAR 013, EAR 008, EAR 022, EAR 041 / HED 021, HED 027 / NKL 001 (mono) Neckloop Charging Kits: BAT KT7 (US), BAT KT8 (Intl)			
	RoHS, WEEE, FCC, IC			
Warranty:	5 years, parts and labor (90 days on accessories)			

Electrical Specifications (1 kHz Sine Wave reference) Typical				
Power	52 mW (max) into 32 ohms (Stereo)			
output:				
ALC:	Automatic Level Control; Limits maximum output level			
Distortion:	1% THD (max)			
Frequency response				
Tone Control at high position:	2.7dB boost at 5 kHz, 6.7dB cut at 400 Hz.			
Tone Control at mid position:	0.6dB cut at 5 kHz, 2.9dB cut at 400 Hz.			
Tone Control				

*Battery switch must be in correct position for battery type being used

at low 2.0dB cut at 5 kHz, 3.0dB boost at 400 Hz.

**Charging requires optinal charging kit BAT KT7 (USA) or BAT KT8 (International).

position:

***When using **mono** earphones, headphones or neckloop, the balance control on the PockeTalker must be turned fully counterclockwise (to the "mono" position) for best volume and battery life.


Note: The PockeTalker is not sold, advertised, or intended for use as a personal hearing aid. Williams Sound recommends that any individual experiencing hearing loss be examined by a physician/ audiologist to determine the cause before using an amplification device such as the PockeTalker.


Acoustic Specifications for Headphones recommended for use with the PockeTalker are published on our website.

Acoustic Specifications: Typical values in dB SPL (ANSI S3.22-1987)

Lightweight stereo headphone, HED 024: (with 6cc coupler)

Test	Tone High	Tone Mid	Tone Low
SSPL90 (Max)	132.2	132.1	131.9
HF Ave (SSPL90)	124.0	124.0	123.9
HF Ave full on gain (60 dB)	45.2	41.7	38.6
Ref test gain	45.2	41.7	38.6

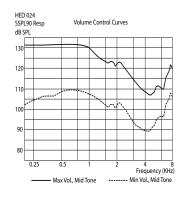


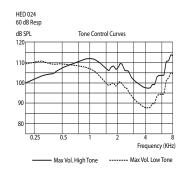

Wide range stereo headphone, HED 036: (with 6cc coupler)


Test	Tone High	Tone Mid	Tone Low
SSPL90 (Max)	128.5	125.7	130.0
HF Ave (SSPL90)	123.4	123.1	124.6
HF Ave full on gain (60 dB)	44.2	40.6	39.3
Ref test gain	44.2	40.6	39.3




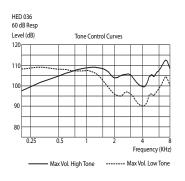

Mini stereo headphone, EAR 042: (with 2cc coupler)

Test	Tone High	Tone Mid	Tone Low
SSPL90 (Max)	124.5	124.5	124.4
HF Ave (SSPL90)	118.9	118.9	118.9
HF Ave full on gain (60 dB)	39.0	35.7	32.8
Ref test gain	39.0	35.7	32.8

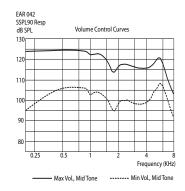

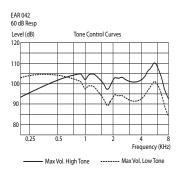

Note: The Pocketalker® 2.0 is not sold, advertised or intended for use as a personal hearing aid. Williams Sound recommends that any individual experiencing hearing loss be examined by a physician/audiologist to determine the cause before using this device.

NOTE: SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE

Domestic Sales

Williams Sound 10300 Valley View Rd Eden Prairie, MN 55344 Ph: 800-328-6190 / 952-943-2252

FAX: 952-943-2174

Email: info@williamssound.com

International Sales

International Sales Department Williams Sound 10300 Valley View Rd Eden Prairie, MN 55344 USA Phone: +1 952 943 2252

Fax: +1 952 943 2174

Email: info-intl@williamssound.com



800.843.3544 / info@williamssound.com / www.williamssound.com